

Canon



Aquilion *Start*

Quality Care for All

Aquilion *Start*



Start with high-end technology

Aquilion Start is packed with premium technology originally developed for our high-end CT systems. It is carefully crafted and adapted to the changing requirements of healthcare organizations, including the need for faster workflows, universal accessibility, and sound economics.

Quality care for management



Opening up new opportunities

With a minimum installation space of only 9.8 m² and minimal power consumption, Aquilion Start fits even into most bucky rooms. Have you ever considered the clinical advantages 3D CT imaging could provide over conventional X-ray?




Quality care for patients

As healthcare becomes more accessible, the number of patients requiring CT scans has significantly increased. Aquilion Start incorporates advanced technologies to optimize patient care and accelerate clinical decision making. Innovative high-end features ensure that you can acquire premium-quality images routinely at very low patient dose.




Integrated dose reduction

Canon Medical's ADR (Adaptive Iterative Dose Reduction) 3D Enhanced is fully integrated into the system's workflow, taking the guesswork out of optimizing patient dose. The exposure dose is automatically reduced by up to 75%.

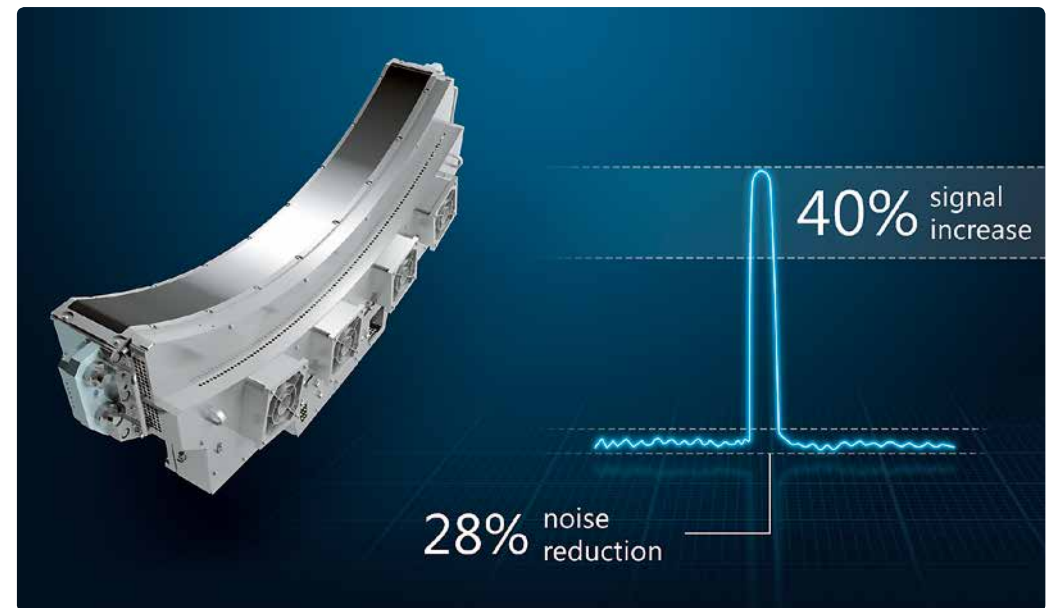

More efficient imaging

PUREVISION Detector technology, which was introduced in our premium systems, has been adapted for Aquilion Start to achieve the perfect balance of image quality and patient dose for outstanding speed in routine clinical practice.

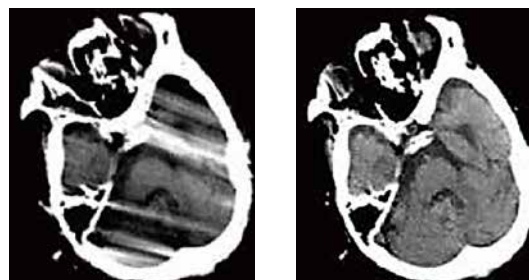


Original

APMC


Advanced motion correction

Canon Medical's APMC (Advanced Patient Motion Correction) technology provides clear, artifact-free images even if the patient moves slightly while being examined.

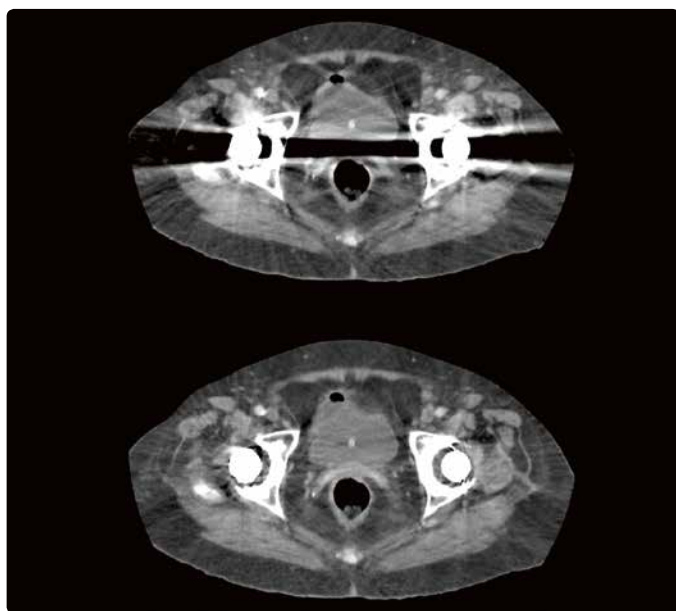


Metal artifact reduction

Dose-neutral SEMAR (Single Energy Metal Artifact Reduction) reduces artifacts caused by metal for better visualization of implants, supporting bone and adjacent soft tissue.

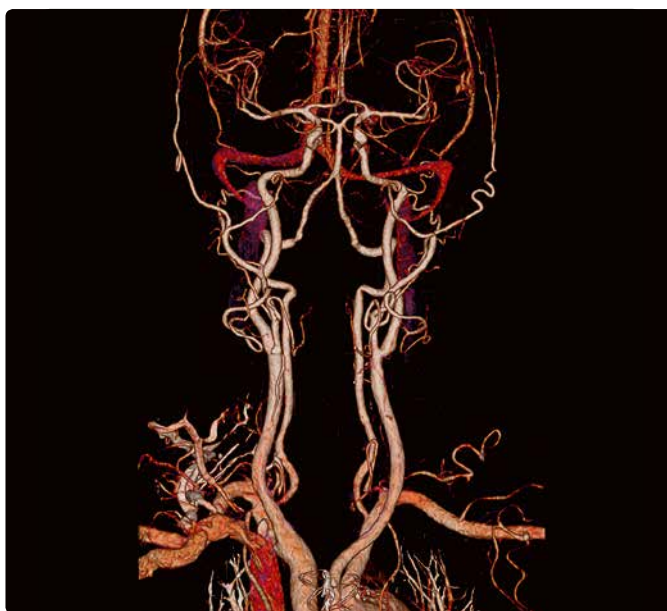

Bone-free imaging*

Subtraction CTA enables superior visualization in CT angiography with true subtraction of bone and calcifications.

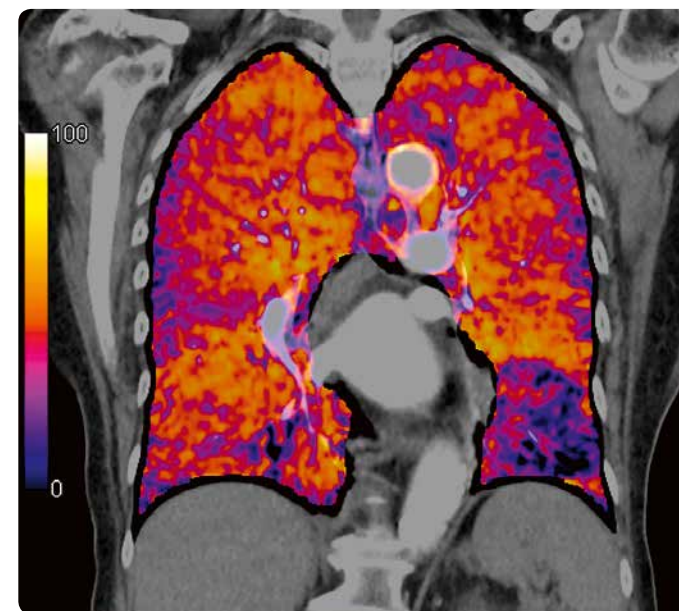


Improved lesion detection*

Iodine Mapping provides clearly defined perfusion with color blood flow maps as a result of advanced registration and subtraction.

*optional

Quality care for operators


Aquilion Start's gantry features design innovations to improve the scanning experience for patients while providing excellent operability and work safety for the technician. The system's wide bore and couch allow comfortable scanning for even the largest of patients. The couch top can be lowered to a minimum height of 31.2 cm for facilitating patient transfer from a wheelchair or stretcher.

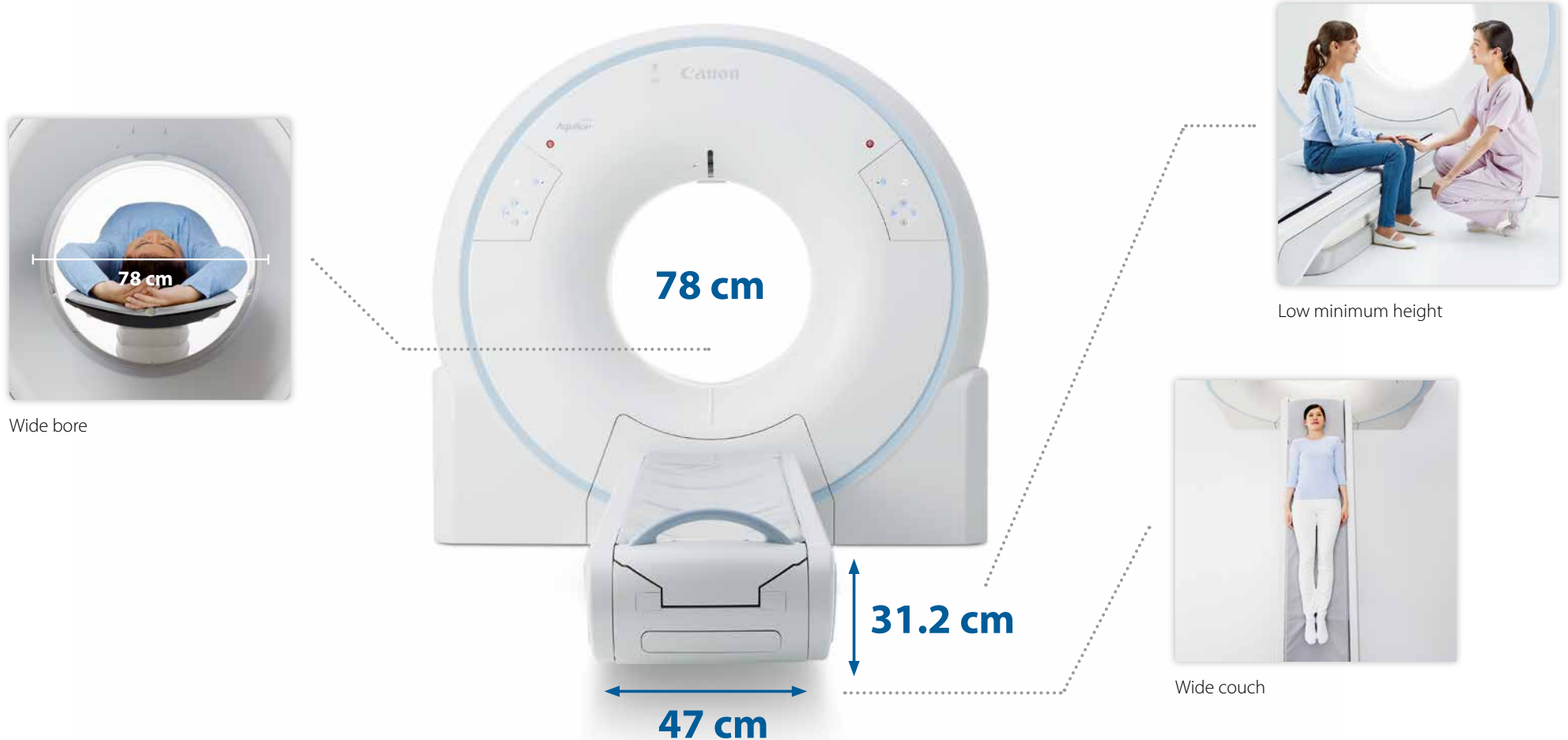



Main specifications		
Detector		PUREVISION detector technology
		16/32 ¹ slices
Gantry	Rotation time	15 s / 1.0 s / 0.75 s
	Bore aperture	78 cm
	Tilt	Flex e-Tilt
Patient couch	Max. load	220 kg
	Max. scan range	183 cm (long couch)
Reconstruction speed		Up to 15 fps
Reconstruction	Iterative reconstruction	AIDR 3D Enhanced
	Power capacity	50 kVA
Installation	Min. installation space	9.8 m ² (short couch) ¹
	Image quality	Spatial resolution

¹ Option ² For reference

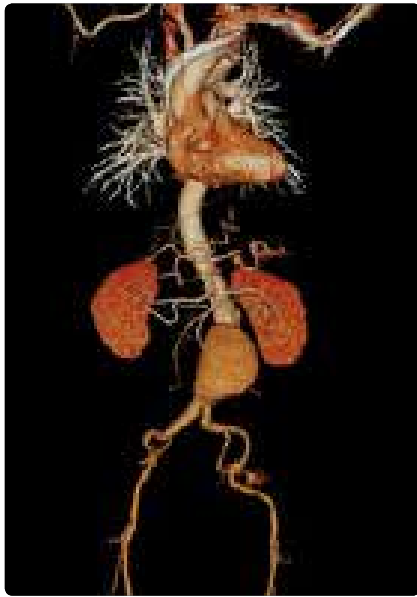

Zero-click workflows

Fully automated, protocol-integrated applications make complex exams easy and provide robust and reproducible results. Zero-click algorithms enable superior image quality without the need for manual interaction, saving valuable time while improving consistency of examinations.



Single-press positioning

With its simplified gantry panels, Aquilion Start streamlines operation and speeds up procedures. Complete patient positioning, for instance, can be done with one single push of a button.



Aquilion Start

Canon

CANON MEDICAL SYSTEMS EUROPE B.V.

<https://eu.medical.canon>

©Canon Medical Systems Corporation 2019. All rights reserved.
Design and specifications subject to change without notice.
Model number: TSX-037A
MCACT0342EC 2019-02 CMSE/Printed in Europe

Canon Medical Systems Corporation meets internationally recognized standards for Quality Management System ISO 9001, ISO 13485. Canon Medical Systems Corporation meets the Environmental Management System standard ISO 14001.

Aquilion, ^{SURE}Exposure, SEMAR and Made for Life are trademarks of Canon Medical Systems Corporation.

Clinical results may vary due to clinical setting, patient presentation and other factors. The clinical images shown in this document are for reference only. They may not have been acquired using Aquilion Start.

Due to local regulatory processes, some of the products included in this brochure may not be commercially available in each country or may require the purchase of additional options. Please contact your sales representative for the most current information.

Made For life